



Trouble-Shooting Guide for ID-Cards

Quick Reference Chart




**Bio-Rad
Laboratories**

Clinical
Diagnostics Group

Website www.bio-rad.com/diagnostics **Australia** +61(2)9914-2800 **Austria** +43-1-877-8901 **Belgium** +32(3)710-53-00 **Brazil** +55(31)3689-6600
Canada +1 514 334-4372 **China** +86-21-61698500 **Czech Republic** +420-241-430-532 **Denmark** +45-4452-1000 **Finland** +358-9-804-22-00
France +33 1 47 95 60 00 **Germany** +49(0)89-318-840 **Greece** +30-210-7774396 **Hong Kong** +852-2789-3300 **Hungary** +36-1-459-6100
India +1-800-180-1224 **Israel** +972-3-9636050 **Italy** +39-02-216091 **Japan** +81-3-6361-7070 **Korea** +82-2-3473-4460 **Mexico** +52(55)5488-7670
The Netherlands +31-318-540666 **New Zealand** +64-9-415-2280 **Norway** +47-23-38-41-30 **Poland** +48-22-3319999 **Portugal** +351-21-472-7700
Russia +7-495-721-14-04 **Singapore** +65-6415-3170 **South Africa** +27-11-442-85-08 **Spain** +34-91-590-5200 **Sweden** +46-8-555-127-00
Switzerland +41 26 674 55 05/06 **Taiwan** +886-2-2578-7189 **Thailand** +662-651-8311 **United Kingdom** +44(0)20-8328-2000



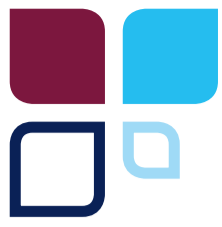

Problems	Possible Causes	Remedy
False/Unexpected Positive Results	▪ ID-Cards showing dry gel	▪ Review storage temperature and conditions
	▪ Expired ID-Card	▪ Discard expired ID-Card and repeat the test with non-expired ID-Card
	▪ Seal damaged	▪ Repeat the test with properly sealed ID-Card
	▪ Air bubbles in the gel	▪ Visual check for bubbles prior to testing
	▪ Direct Antiglobulin Test (DAT) positive cells	▪ Repeat using monoclonal ID-Card or warm-washed cells
	▪ Polyagglutinable cells	▪ Repeat using monoclonal ID-Cards
	▪ Interference by enzyme autoagglutinins	▪ Repeat using monoclonal ID-Cards
	▪ Carry-over of antisera between microtubes	▪ Check pipetting technique ▪ Pipette cells in reverse order ▪ Ensure ID-Card is upright when opened ▪ If droplets under the foil, centrifuge the ID-Card prior to use
	▪ Incorrect cell suspension	▪ Repeat with correct cell suspension
	▪ Wharton's Jelly - cord blood	▪ Wash cells and repeat
Equivocal Results	▪ Centrifugation problems (ID-Card holder stuck, speed/time issues)	▪ Report to your field service contact point
False/Unexpected Weak or Negative Results	▪ Plasma/serum omitted from test procedure	▪ Repeat test - check level visually for plasma/serum addition
	▪ Wrong volume(s) used	▪ Check box insert for correct procedure
	▪ Deterioration of antibody (poor storage, natural phenomenon)	▪ Fresh samples should be used for optimal results
	▪ Weak/missing isoagglutinins	▪ Review clinical details (age of patient, medical diagnosis, therapy, etc.)
	▪ Weak/missing antigens	▪ Review clinical details (age of patient, medical diagnosis, therapy, etc.) ▪ Subgroup determination (ABO) ▪ Test for weak D/partial D
Double Cell Population Mixed Field Reaction	▪ Contamination of reagent cells	▪ Repeat with fresh set of cells ▪ Review laboratory procedure
	▪ Post-transfusion	▪ Review clinical history/treatment
	▪ Post-Stem Cell/Bone Marrow Transplant	▪ Review clinical history/treatment
	▪ Large fetomaternal hemorrhage	▪ Appropriate further investigation
	▪ Subgroups (A ₃ , A _{weak} , etc.)	▪ Appropriate further investigation
	▪ Chimerism	▪ Appropriate further investigation
	▪ Sample contamination	▪ Repeat with fresh sample (review phlebotomy/laboratory procedures)
	▪ Fibrin residues (trapping red blood cells)	▪ Serum/plasma should be well centrifuged before use
	▪ Weakened antigens	▪ Review clinical details e.g. in association with leukemia, old age, newborns etc.

Notes

- For specific problems in ABO or Rh typing, refer to the appropriate Quick Reference Chart for guidance.
- All laboratory procedures should be subject to internal quality control as part of Good Laboratory Practice.
- Report persistent problems to your local customer product support team, who will take the appropriate action required to resolve any difficulties.
- This chart is not necessarily comprehensive.